

Cells on micropatterned substrates

Glass coverslip  clean & activate w/UV

PLL-PEG  parafilm

PLL-PEG coated covers. 

a) incubate $\frac{1}{2}$ -1h

b) clean w/pure water or PBS
dry w/ N_2
can be stored!

Quartz mask 

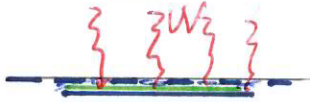
clean & activate w/UV

coverslip on mask 

a) small volume water $V = 5 \mu m \cdot A$

b) place coverslip, capillary force pulls it close

c) place under lamp & illuminate
can be stored!

Illuminate 

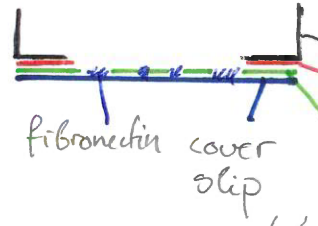
PLL-PEG-
Patterned coverslip 

Fibronectin 

a) incubate 1h in 25 $\mu g/ml$ fibronectin in 100 mM $NaHCO_3$

b) rinse w/PBS & dry

Willco Well

 WW petri dish
double-sided adhesive
PLL-PEG
fibronectin cover slip

a) Glue cover slip to Willco well petri dish.

seed cells



a) pipette suspended cells into dish
b) add media

c) place in microscope in incubator

Image



Analyse

weeks

1h

1day